The advertisement features a muscular man and woman in the background. The man is shirtless, showing his abdominal muscles, and the woman is wearing a black sports bra and a black waist trainer with an orange band. In the foreground, a black abdominal exerciser with orange patterns is shown. The exerciser has a central circular dial and eight segments with orange lines. The text 'GYM form TOTAL ABS PLATINUM SPORT' is prominently displayed in the upper right. At the bottom, there are four language options for the instruction manual, each with a small flag icon.

GYM[®] form

TOTAL ABS PLATINUM

SPORT

Instruction Manual

Manual de instrucciones

Mode d'emploi

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






International patents pending. All imitations will be prosecuted.
Gymform® Total ABS Platinum PLATINUM is a registered EU/CTM trade mark.

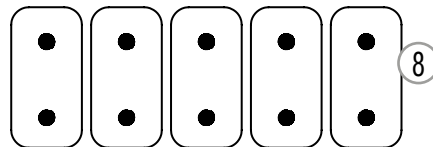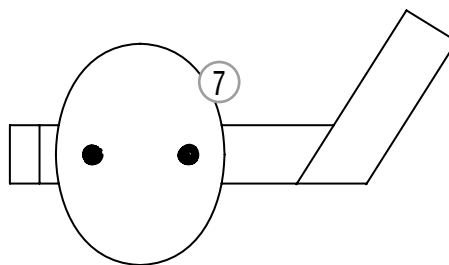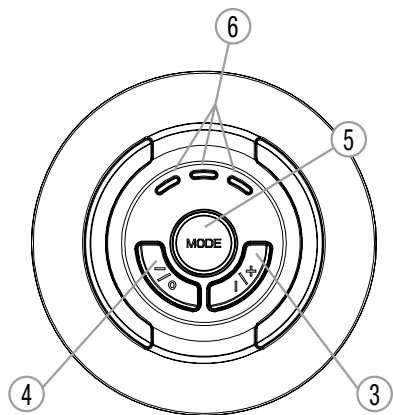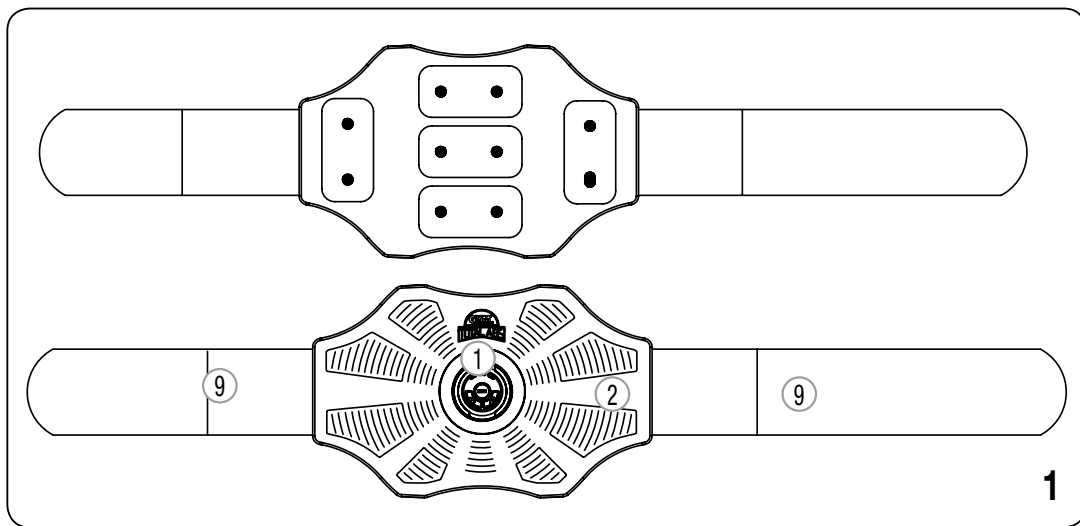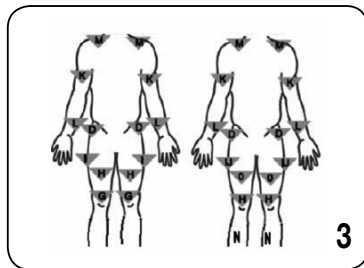
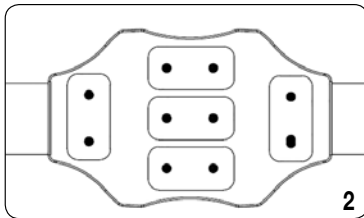
Programme 1 – MASSAGE
Programma 1 – MASAJE
Programme 1 – MASSAGE
Programm 1 – MASSAGE
Programma 1 – MASSAGGIO
Programma 1 – MASSAGEM
Programma 1 – MASSAGE

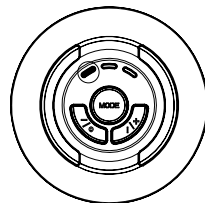
Programme 4 – POWER SCULPT
Programma 4 – ESCULPIR MÚSCULO
Programme 4 – SCULPTER LES MUSCLES
Programm 4 – MUSKEL AUSBILDEN
Programma 4 – SCOLPIRE I MUSCOLI
Programma 4 – ESCULPIR MÚSCULOS
Programma 4 – VORMGEVEN AAN DE SPIER

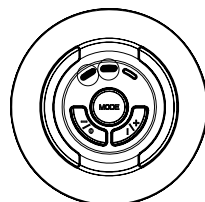
Programme 2 – FAT BURN
Programma 2 – QUEMAGRASAS
Programme 2 – BRÛLEUR DE GRAISSES
Programm 2 – FETTVERBRENNUNG
Programma 2 – BRUCIAGRASSI
Programma 2 – QUEIMA-GORDURA
Programma 2 – VETTERBRANDUNG

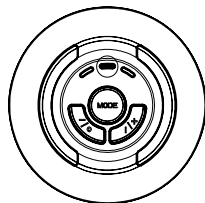
Programme 5 – TOTAL EXERCISE
Programma 5 – ENTRENAMIENTO COMPLETO
Programme 5 – ENTRAÎNEMENT COMPLET
Programm 5 – KOMPLETTES TRAINING
Programma 5 – ALLENAMENTO COMPLETO
Programma 5 – TREINO COMPLETO
Programma 5 – VOLLEDIGE TRAINING

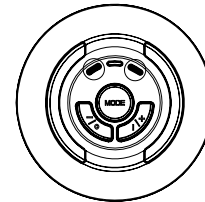Programme 3 – MUSCLE BUILD
Programma 3 – MUSCULACIÓN
Programme 3 – MUSCULATION
Programma 3 – MUSKELAUFBAU
Programma 3 – MUSCOLAZIONE
Programma 3 – MUSCULAÇÃO
Programma 3 – SPIEROPBOUW

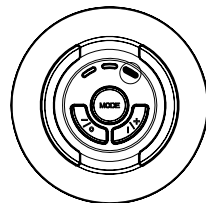
Programme 6 – CELLULITE DRAIN
Programma 6 – DRENAJE CELULÍTICO
Programme 6 – DRAINAGE DE LA CELLULITE
Programma 6 – CELLULITE-DRAINAGE
Programma 6 – DRENAGGIO ANTICELLULITE
Programma 6 – DRENAGEM CELÚLTICA
Programma 6 – ANTI-CELLULITISDRAINAGE

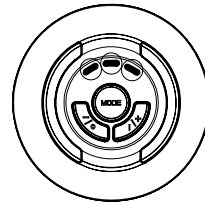
Programme 7 – STRENGTH TRAINING
Programma 7 – ENTRENAMIENTO DE FUERZA
Programme 7 – ENTRAÎNEMENT DE FORCE
Programm 7 – KRAFTTRAINING
Programma 7 – ALLENAMENTO DI FORZA
Programma 7 – TREINO DE FORÇA
Programma 7 – POWER SCULPT
(KRÄCHTTRAINING)

Flashing light
Luz intermitente
Lumière clignotante
Blinklicht
Luce intermittente
Luz intermitente
Knipperlicht

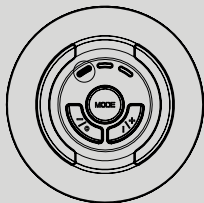
Programme 8 – FIRMING AND TONING
Programma 8 – REAFIRMANTE Y TONIFICANTE
Programme 8 – RAFFERMISSEMENT ET TONIFICATION
Programm 8 – STRAFFEND UND STÄRKEND
Programma 8 – RASSODANTE E TONIFICANTE
Programma 8 – REFIRMANTE E TONIFICANTE
Programma 8 – VERSTEVIGEND EN SPIERVERSTERKEND

Flashing light
Luz intermitente
Lumière clignotante
Blinklicht
Luce intermittente
Luz intermitente
Knipperlicht

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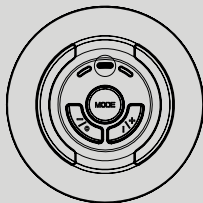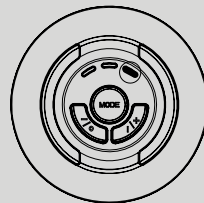Congratulations on your purchase of the **Gymform® Total ABS Platinum Platinum**. Your safety is our primary concern, so please ensure that you familiarize yourself with all the information in this manual before using your **Gymform® Total ABS Platinum**. Anyone who intends to use this device must read and understand the instructions provided. An informed user is a safe user. Thank you for purchasing the **Gymform® Total ABS Platinum**. We hope you enjoy your sessions.

WARNINGS

Before embarking on any exercise programme, including the **Gymform® Total ABS Platinum**, consult your general practitioner and your physiotherapist. Ensure that you read the instructions, warnings and precautions carefully before using your **Gymform® Total ABS Platinum**.

Gymform® Total ABS Platinum is designed for healthy adults. However, the device should be used with care. Read the following guidelines and instruction manual before using the product.

- Applying electrodes in the thorax area may increase the risk of cardiac fibrillation.
- Using the **Gymform® Total ABS Platinum** while connected to a high-frequency surgical unit may cause burns at the contact points of the **Gymform® Total ABS Platinum** and may result in damage to your

Gymform® Total ABS Platinum.

- Operation of a short-wave therapeutic device or microwave oven in the general vicinity (e.g. 1 metre away) may cause instability in the output power of the **Gymform® Total ABS Platinum**.

- This is a powerful electrical muscle stimulation device and may cause injury if you start with the most advanced exercise programme right away. Start with a lower level and gradually increase the intensity.

- Do not use on children or animals.

- Do not use while driving or sleeping.

- Do not use in the bathroom or in other places with a high amount of moisture or humidity.

- The batteries should be removed from the device before it is disposed of.

- The batteries should be safely disposed of.

- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe

way and understand the hazards involved

- Children shall not play with the appliance

- Cleaning and user maintenance shall not be made by children without supervision

- Non rechargeable batteries are not to be recharged

- Exhausted batteries are to be removed from the appliance and safely disposed of

- If the appliance is to be stored unused for a long period, the batteries are removed the supply terminals are not to be short circuited

DO NOT use Gymform® Total ABS Platinum:

- If you have an electronic implant (such as a pacemaker) or suffer from another heart problem.

- If you are, or think you might be, pregnant.

- If you have cancer, epilepsy or are under medical supervision for a cognitive disorder.

- If you are breastfeeding.

- If you have a high fever.

Women, before using Gymform® Total ABS Platinum, you should wait:

- At least six months after the birth of your child (after consulting your doctor).

- A month after having an intrau-

terine device (IUD) fitted.

- At least three months after having a caesarean section (after consulting your doctor).
- Until after the first days of your period.

Before using Gymform® Total ABS Platinum, you should obtain permission from your doctor or physiotherapist:

- If you suffer from back problems. In this case, the device should only be used at low intensity levels.
- If you suffer from a serious disease or injury not mentioned in this section.
- If you were operated on recently.
- If you are diabetic and take insulin.

When putting on the belt, always remember:

- To place the large belt ONLY on your abdomen, obliques and lower back.
- To place the small belt ONLY on the areas indicated below.
- Not to place the belt on your neck, on your chest or on your back in such a way that electricity passes through your heart (e.g. an electrode on your chest and another on your back), genitals, scars, moles, varicose veins or head.
- Not to place the belt on recent

scars, inflamed skin, infected areas or areas that are susceptible to acne, thrombosis or other vascular problems (e.g. varicose veins), or any part of your body with limited sensitivity.

- Do not place on areas with injuries or limited mobility (e.g. fractures or sprains).
- Not to place the belt directly on metal implants. Place on the closest muscle.

Possible Adverse Reactions

There have been isolated reports of skin reactions caused by the use of electro-stimulation devices, including allergies, prolonged skin redness and acne. If any of these symptoms occur, switch off the device and consult your doctor.

If you have never exercised before or have not exercised for a long time:

Consult your general practitioner.

Precautions:

Do not use **Gymform® Total ABS Platinum** if you are drowsy or while sleeping.

Do not use **Gymform® Total ABS Platinum** if your body is in contact with water.

- Keep out of reach of children.
- Do not use **Gymform® Total ABS Platinum** in combination with any other device that trans-

mits electrical current to your body or with any other machine.

- Do not connect the metal snap fasteners or electrodes to any other object.
- Do not modify the device.
- Use only the electrodes supplied by Industex.
- Do not attach two electrodes together.
- If any of the electrodes is not attached properly, you will notice an itching sensation.

Note:

- The product may malfunction as a result of electrostatic disturbances. If this occurs, simply press the RESET button (16) to resume normal operation.

COMPONENTS:

- 1 Detachable **Gymform® Total ABS Platinum** unit with 9 programmes.
- 1 **Gymform® Total ABS Platinum** belt
- 1 Special belt for working the leg and arm muscles
- 5 Self-adhesive electrodes
- 1 Instruction manual.
- 1 Table of fitness targets
- 1 Sports bag
- 1 Diet booklet

PRODUCT DESCRIPTION (Fig. 1):

1. Gymform® Total ABS Platinum unit
2. Gymform® Total ABS Platinum belt
3. + / I button (to turn device on and increase intensity level)
4. - / o button (to turn device off and reduce intensity level)
5. Mode button (programme selection)
6. Programme indicator lights
7. Arm and leg belt
8. Electrodes

9. Adjustable fastener

Putting on the Belt:

Set-up:

Open the battery compartment cover (batteries not included) on the back of the unit (1) and insert three AAA batteries. To remove the battery compartment cover, insert the screwdriver into each screw head and turn until the screw is released. After inserting the batteries, replace the cover and position the screws. Then reinsert the screwdriver into the screw head of each one and turn until it fits tightly. Insert the three batteries in accordance with the polarity marks. The polarity is indicated on the battery compartment for proper placement.

Putting on the Belt:

Attach the Gymform® Total ABS Platinum (1) to the front of the leg and arm belt (2) using the two snap fasteners. Ensure that they are both well fastened or the belt will not work properly.

Attach the five electrodes (8) to the inside of the belt (2) as shown in Figure 2. Use the snap fasteners so that the electrodes are firmly attached.

Remove the protective plastic backing from the belt electrodes. Do not throw the protective plastic backing away, as you will need it later.

Place the belt around your abdomen so that one electrode is on your upper abs, one is on your middle abs, one is on your lower abs, and two are on your obliques. Ensure that all the parts of the belt fit properly and the electrodes are in contact with the skin. Keep the belt firmly in place using the adjustable fastener.

When you have finished using the belt, to ensure the electrodes do not stick to the skin, remove them one at a time by pulling them from the edge and then replace the protective plastic backing on the electrodes. This will help extend the useful life of the electrodes (8).

Programme Selection

Press the + / I button (3) to turn on the unit (1).

When **Gymform® Total ABS Platinum** is first switched on, it starts at Programme 1 with the intensity at 0. Use the MODE button (5) to select other programmes. Press the MODE button (5) once to change to Programme 2, once again for Programme 3, and so on. After reaching Programme 6, it will return to Programme 1. The programme selected is clearly shown by the programme indicator lights (6) on the **Gymform® Total ABS Platinum** unit. See the table on the second page to identify the programmes.

Adjusting the Intensity Level

To reduce the intensity level, briefly and repeatedly press the - / o button (4). To increase the intensity level, briefly and repeatedly press the + / I button (3). There are 10 intensity levels to choose from. If this is your first electro-stimulation session, we recommend that you start at Level 1 for the first 10 minutes and then gradually increase the intensity as you become more familiar with the device.

If you press - / o (4) for two seconds, the device will turn off.

Automatic Power Off

For your safety, the automatic period for each programme is 10 minutes. **Gymform® Total ABS Platinum** will automatically switch itself off 10 minutes after being switched on, regardless of the level or programme selected.

If you wish to finish your session before the end of the programme time, press the - / o button (4) for one second until the device switches off.

Automatic Reset

When you finish the programme, the device will switch itself off automatically. When you switch it on again, you will have to programme it again.

Using the Leg and Arm Belt

Attach the Gymform® Total ABS Platinum unit (1) to the front of the leg and arm belt (7) using the two snap fasteners. Ensure that they are firmly attached.

Attach one electrode (8) to the belt using the snap.

Remove the protective plastic backing from the electrode. (Make sure you replace it after use.)

Position the belt on the area you wish to work on (see figure below) and ensure that the belt is well secured with the quick opening and closing system so that the electrodes are in contact with the skin.

To choose the programme and intensity level, see the section above.

AREAS WHERE THE LEG AND ARM BELT MAY BE USED (Fig. 4)

Front

- G) Lower front thigh
- H) Upper front thigh
- I) Outer upper thigh (front)

Rear

- J) Outer upper thigh (back)
- K) Arms: biceps and triceps
- L) Forearms
- M) Trapezius muscles
- N) Calves
- O) Inner thighs

WE RECOMMEND THAT YOU REPLACE THE ELECTRODES AFTER EVERY 30 USES. IF YOU NOTICE THE ELECTRODES HAVE PREMATURELY LOST THEIR ADHESIVE POWER, WET THEM WITH A DAMP SPONGE.

NOTE: If you feel heat in the areas where you have attached the electrodes or if the contractions become weak during the session, it is likely that the batteries need to be replaced. Remove them and insert new ones. If this does not solve the problem, the electrodes have probably reached the end of their useful life. Switch off the **Gymform® Total ABS Platinum**, take it off, fit new electrodes to the belt

and then resume your session. For your convenience, **Gymform® Total ABS Platinum** will automatically set the intensity level to 0 if you change the programme.

TECHNICAL SPECIFICATIONS:

Temperature: 0~+35°C

<85% RH

700-1060 hPa

Continuously operating device.

Power supply: 3 1.5V AAA batteries (not included)

Voltage: 3 x 1.5V (4.5V total)

Program	Pulse with (us)	Frecuency (Hz)	Voltage (V)
1	200	70	38 V
2	200	60	38 V
3	200	50	38 V
4	200	50	38 V
5	200	50	38 V
6	200	50	38 V
7	122	10	38 V
8	122	50	38 V
9	108	2	38 V

Values measured at 4 kΩ

Maximum current of the electrodes supplied:

Dimensions: Φ 100 × 23.7 mm

Weight: 61 g (main unit without batteries)

Operating voltage: DC 4.5V (AAA x 3)

Safety classification based on the type of protection against electric shock: internally powered device with a type BF applied part.

ISO 7000-1641 User Guide

CLEANING AND MAINTENANCE

DO NOT WASH THE **Gymform® Total ABS Platinum** UNIT OR THE BELT: WIPE WITH A DRY OR DAMP CLOTH.

Remove the battery from your Gymform® Total ABS Platinum if you do not plan to use it for an extended period of time.

MAINTENANCE AND CLEANING OF ADHESIVE ELECTRODES

- Clean by wiping with wet fingers.
- Do not immerse in water.
- Keep in a cool, dry place at room temperature.
- Do not remove the insulating label.
- After using the device, always replace the protective plastic backings on the adhesive electrodes.
- Store in the bag supplied.

PRECAUTIONS WITH BATTERIES

Ensure the batteries are properly installed before turning the unit on.

Turn off the unit to change the batteries.

Remove the batteries if you are not going to use the device for an extended period of time.

Take the necessary precautions when disposing of the batteries. Do not expose them to fire.

DISPOSAL OF MATERIALS

A symbol of a crossed-out wheeled bin means you should find out about and follow local regulations about disposing of this kind of product.

Do not dispose of this product as you would other household waste. Dispose of this appliance in accordance with the corresponding local regulations.

Electrical and electronic appliances contain hazardous substances that can have harmful effects on the environment and/or human health and should be recycled properly.

The battery should be disposed of separately. To protect the environment, all batteries, including rechargeable ones, should not be disposed of as you would other household waste. Take them to a disposal centre and ensure it complies with applicable legislation on battery disposal.

Quality Warranty:

This product is covered by a warranty against manufacturing defects subject to the time periods stipulated by the legislation in force in each country.

This warranty does not cover damages resulting from inadequate use, negligent commercial use, abnormal wear and tear, accidents or improper handling.

Australia & New Zealand only

This product has the benefit of certain Consumer Guarantees. These are prescribed by the Australian Consumer Law & the New Zealand Consumer Guarantees Act 1993, both of which provide protection for consumers. There is no express warranty for this product in Australia or New Zealand. The above paragraph refers to other countries.

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ISL GmbH
Best Direct (International) Benelux B.V.
Best Direct International Ltd,
ISL ITALY S.r.l.